

Dansylated Monoamines

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Methanol

Eluent B: Water mixture

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| | 1 | 0 | 55 | 45 |
| | 2 | 20 | 90 | 10 |

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: Fluorimetric @ nm (ambient)

Detector Info: The following Fluorometric filters were employed: 7-60/7-54 and 3-72/4-76

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

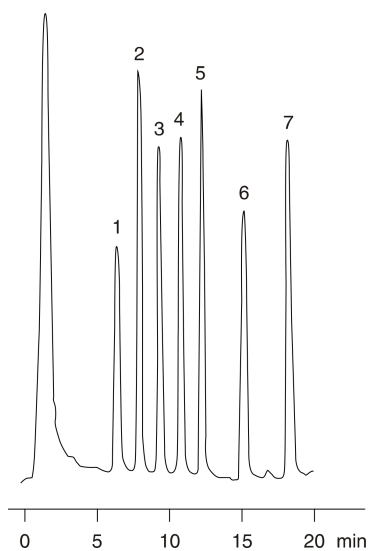

ANALYTES:

- 1 Ethylamine
- 2 Propylamine
- 3 Butylamine
- 4 Pentylamine
- 5 Hexylamine
- 6 Octylamine
- 7 Decylamine

