

## Tocopherols - SphereClone vs. Spherisorb Comparison

**Column:** SphereClone<sup>™</sup> 5  $\mu$ m ODS(2) 80 A, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4144-E0

**Elution Type:** Isocratic

**Eluent A:** Acetonitrile

**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 295 nm (ambient)

**Analyst Note:** SecurityGuard<sup>™</sup> Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

**SphereClone<sup>™</sup>**

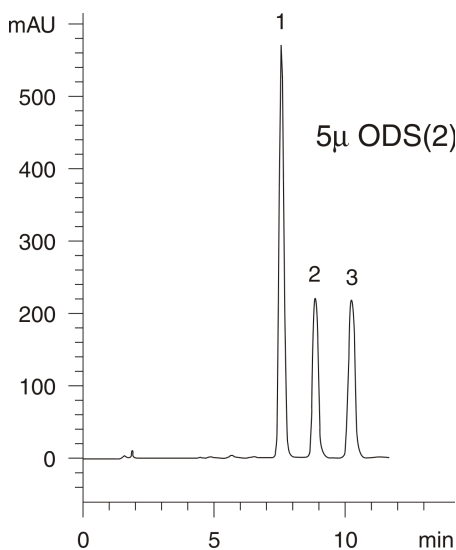

### ANALYTES:

- 1 delta-Tocopherol
- 2 gamma-Tocopherol
- 3 alpha-Tocopherol

