

## Glucocorticoids - SphereClone vs. Spherisorb Comparison

**Column:** SphereClone™ 5 μm ODS(1) 80 A, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4143-E0

**Elution Type:** Isocratic

**Eluent A:** Acetonitrile

**Eluent B:** Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

**SphereClone™**

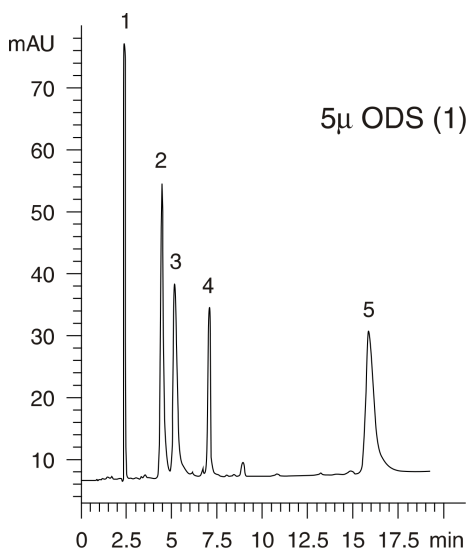

### ANALYTES:

- 1 Uracil
- 2 Prednisolone
- 3 Betamethasone
- 4 Prednisolone 21-acetate
- 5 Betamethasone 17-valerate

