

## Sugars in Honey on Rezex ROA

**Column:** Rezex™ ROA-Organic Acid H+ (8%), LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-0138-K0

**Elution Type:** Isocratic

**Eluent A:** 0.01 N H<sub>2</sub>SO<sub>4</sub>

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index (RI) @ nm (ambient)

**Analyst Note:** Standard solution (A) vs. Honey sample (B)

SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-H 4 x 3.0mm, 10/Pk Part No.: AJ0-4490


- Holder Part No.: KJ0-4282

**REZEX™**  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX™**

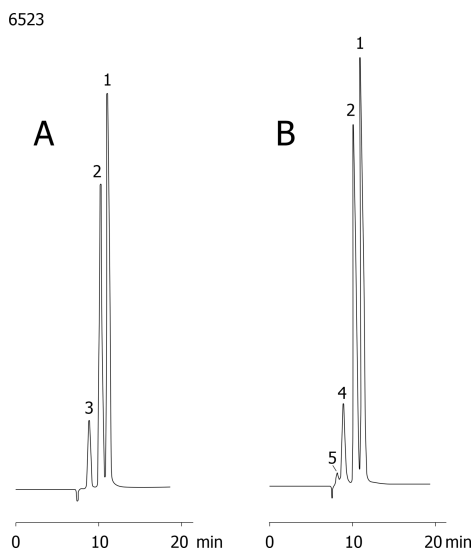

### ANALYTES:

- 1 Fructose
- 2 Glucose
- 3 Sucrose
- 4 Disaccharide
- 5 Trisaccharide

