

Drugs of Abuse

Column: Zebron™ ZB-1, GC Cap. Column 15 m x 0.25 mm x 0.25 µm, Ea

Phase: 100% Dimethylpolysiloxane

Dimensions: 15 meters x 0.25 mm x 0.25 µm

Order No: 7EG-G001-11

Oven Profile: 70°C to 200°C @ 35°C/min, to 250°C @ 23°C/min, to 325°C @ 20°C/min


Carrier Gas: Constant Pressure Helium, 1 mL/min

Injection: Split :1 1 µL @ 250°C

Detection: Mass Spectrometer (MS) (280°C)



Products used in this application:

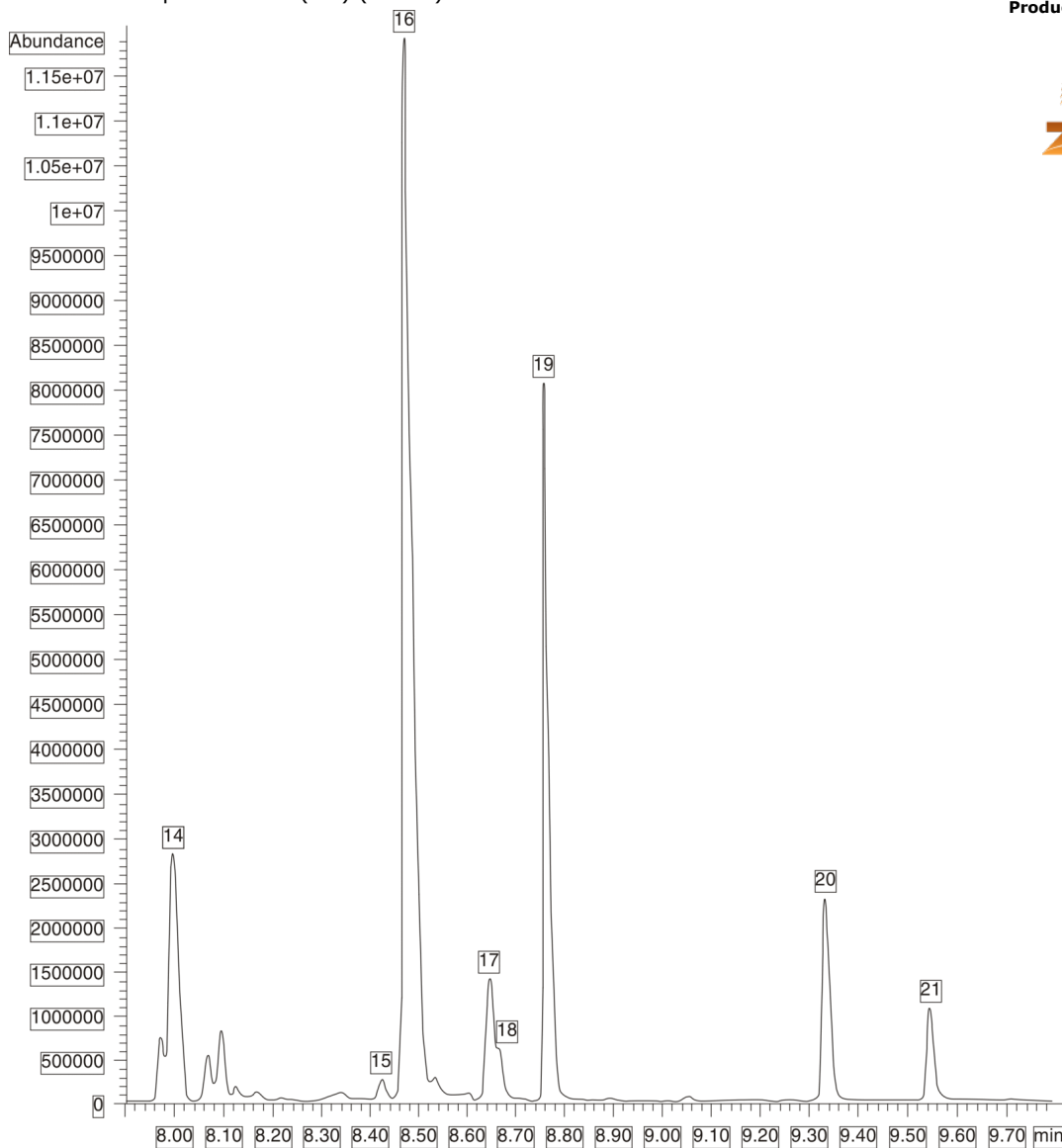

Drugs of Abuse

ANALYTES:

- 1 Andyroecgonine methyl ester
- 2 Methamphetamine
- 3 Amphetamine
- 4 Ecgonine methyl ester
- 5 Ecgonine ethyl ester
- 6 Clonidine
- 7 Methadone
- 8 Cocaine
- 9 Cocaethylene
- 10 Norcocaine
- 11 Norcocaethylene
- 12 Codeine
- 13 Benzoylecgonine
- 14 Norbenzoylecgonine
- 15 Morphine
- 16 m-Hydroxycocaine
- 17 6-Acetylmorphine
- 18 Normorphine
- 19 p-Hydroxycocaine
- 20 m-Hydroxybenzoylecgonine
- 21 p-Hydroxybenzoylecgonine

